PEDANA VIBRANTE

IT

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VIBRATING PLATFORM

GB

USER MANUAL

ATTENZIONE!

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WARNING!

We thank you for the trust placed in us by purchasing a product branded

The contents of this manual are related to the use and maintenance of this equipment, so please keep it for later consultation.

INSTRUCTION MANUAL

Please read carefully before use and keep for future reference. If a product malfunction occurs and its resolution is not indicated in this manual, DO NOT try to fix it on your own but please contact qualified technical assistance. Due to the constant development of product functions and designs, we reserve the right to make changes to the product without prior notice.

Importer and distributor: RICA TRADE S.R.L.

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DECLARATION OF CONFORMITY

The device described in this documentation complies with the legal provisions regarding safety at the time of production.

MANUFACTURER'S DATA

Company name	Rica Trade s.r.l.
Company address:	LUNGOTEVERE DEI MELLINI, 44, SCALA VALADIER 2/B 00193 ROMA
VAT number	09785201006
E-mail	infoprodotti@ricatrade.it
Website	www.risparmiocasa.com

PRESENTATION OF THE MANUAL

ON RECEIPT OF THE DEVICE, BEFORE PERFORMING ANY OPERATION, READ THIS MANUAL CAREFULLY.

This manual contains instructions for use and maintenance of the device purchased. The manual consists of several sections, each of which deals with a number of topics.

This manual is intended for the user in charge of using the device and relates to the technical life after its production and sale. In the event that it is subsequently transferred to third parties for any reason (sale, loan for use, or any other reason), the device must be delivered complete with all documentation.

Before starting any operation it is necessary to have at least read the entire manual and then deepened the topic relating to the operations to be carried out. This manual contains proprietary information and may not even be partially provided to third parties for any use and in any form. We declare that the information contained in this manual is consistent with the technical and safety specifications of the device to which the manual refers. A certified copy of this manual is deposited in the technical file of the device, kept at the manufacturer, which does not recognize any documentation that has not been produced, issued or distributed by itself or by its authorized representative. This manual, like all the technical file, will be kept by the manufacturer for the period required by law. After this period, it will be the responsibility and obligation of those who manage the product to ensure that both the product and the documentation comply with the laws in force at the time of the inspection.

WARRANTY

The warranty rules are valid only if the device is used in the conditions of intended use.

Any repair or modification made to the device by the user or by unauthorized companies will invalidate the warranty.

The warranty does not extend to damage caused by inexperience or negligence in the use of the device, or by bad or omitted maintenance.

The products we sell are covered by a warranty with respect to the device under the following conditions:

The warranty is valid for a period of 12 or 24 months depending on the type of buyer (legal or private). The manufacturing company undertakes to replace at its discretion the malfunctioning or incorrectly manufactured parts, only after careful control and verification of bad assembly. The transport and / or shipping costs in the event of incorrect use of the warranty terms are always borne by the buyer. During the warranty period, the replaced products become the property of the manufacturer. Only the original purchaser who has complied with the instructions for normal maintenance contained in the manual can benefit from this guarantee. Warranty liability expires when the original owner relinquishes ownership of the product, or changes have been made to the product. The warranty does not include damage resulting from excessive stress, such as the use of the product after finding an anomaly, from the use of unsuitable operating methods as well as from failure to observe the instructions for use and maintenance. The manufacturer assumes no responsibility for any difficulties that may arise in resale or use abroad due to the provisions in force in the country in which the product was sold. The product or part of the defective product must be delivered to the manufacturer at its own expense for replacement and will not be replaced without proof of purchase (invoice, copy of payment, receipt); otherwise the replaced part will be charged to the buyer.

GENERAL SAFETY RULES

Warnings

If part of the documentation is even partially missing or illegible, consult the manufacturer before carrying out any further operation on the device.

This chapter describes the general safety rules to be observed during any operation carried out with the device. The intervention procedures, described in the following chapters, must be carried out in compliance with both the methods of execution indicated and the general safety rules of this chapter.

The safety rules and the configuration, use and maintenance procedures indicated in this document can also be a complement to the general workplace safety rules that must be respected.

Different nations may have different safety regulations. It is therefore specified that in all cases where the documentation standards are in conflict or reductive with respect to the standards of the country in which the equipment is used, the standards of the country will still have priority over those of the documentation.

THE MANUFACTURER CANNOT BE HELD IN ANY CASE LIABLE FOR

ACCIDENTS OR DAMAGES RESULTING FROM THE INAPPROPRIATE USE OF THE DEVICE, AS WELL AS EVEN PARTIAL NON-COMPLIANCE WITH THE SAFETY RULES AND INTERVENTION PROCEDURES DESCRIBED IN THE DOCUMENTATION.

Failure to comply with the rules of use and the methods of intervention, configuration, use and maintenance of the device contained in the manual also results in the cancellation of the warranty terms.

PRODUCT DESCRIPTION:

The vibrating platform provides a very complete training method, which affects not only muscle mass, but also other aspects related to good physical shape, for example by stimulating blood circulation. The training on the vibrating platform has a significant positive impact on the body, optimally stimulates the metabolism, with beneficial effects for weight loss.

Component list:

Vibrating platform

1. Control Panel
2. Elastic band hook
3. Power cord socket
4. Main switch
5. Support feet
6. Support surface


Accessories:

2 x elastic bands


Control panel:

7. TIME+: to increase the training time
8. Timer and program display: remaining time in minutes, selected program
9. Remote control receiver
10. Speed display: vibration level
11. SPEED+: to increase the vibration level
12. SPEED-: to decrease vibration level
13. AUTO: to select a predefined program
14. MANU: manual program
15. STOP
16. START
17. TIME-: to decrease the time

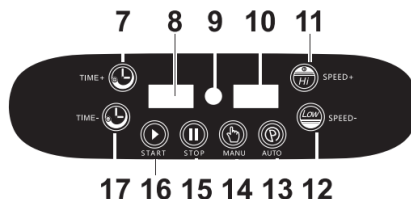Remote controller:

- 18. ►: START
- 19. TIME+
- 20. SPEED-
- 21. M: manual program
- 22. Quick select vibration level 30
- 23. Quick select vibration level 60
- 24. Quick select vibration level 90
- 25. ■: STOP
- 26. AUTO
- 27. TIME-
- 28. SPEED+
- 29. ⏻: STAND BY/OFF
- 30. Remote controller light indicator

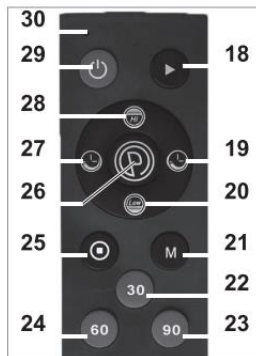TECHNICAL DATA

DIMENSIONS: 59.7X34X11.7CM

WEIGHT: 10.3kg

MATERIAL: Plastic and TPR

POWER SUPPLY: AC220-240V 50-60Hz

POWER (W): 200W

SAFETY CLASS: class I

Maximum capacity: 100kg

IP: IPX0

Noise: 70dB

99 levels of vibration

10 predefined programs

With remote control

Accessories: 2 elastic bands

THE DATA CONCERNS THE DEVICE IN STANDARD CONFIGURATION. THE TECHNICAL DATA REPORTED ARE INDICATIVE. THE MANUFACTURER RESERVES THE RIGHT TO MODIFY THEM APPROPRIATELY WITHOUT NOTICE.

ATTENTION!

VERIFY THAT THE ELECTRICITY NETWORK TO WHICH THE SYSTEM WILL BE CONNECTED COMPLIES WITH THE ESSENTIAL SAFETY REQUIREMENTS IN FORCE IN THE TERRITORY AND THERE IS A LIFE-SAVING SYSTEM OR RCD whose nominal operating differential current does not exceed 30mA.

THE CONNECTION CABLES SHOULD NEVER COME IN CONTACT WITH WATER OR SPLASH. THEY SHOULD NOT BE POSITIONED IN A WAY TO SUFFER DAMAGE. IF THE ELECTRIC CABLES ARE DAMAGED, YOU MUST STOP USING THE DEVICE. THE CABLES MUST NOT BE INTRACTING FOR PEOPLE AND ANIMALS.

THE CONNECTION CABLES MUST NOT BE BURIED, WALLED OR HIDDEN BEHIND CURTAINS OR UNDER RUGS. THEY MUST REMAIN INSPECTIONABLE FOR THE WHOLE TIME OF USE.

WHEN THE DEVICE IS DISCONNECTED FROM THE ELECTRICITY NETWORK, THE CABLE MUST NOT BE ACCESSIBLE, IT COULD BE INVOLUNTARILY PULLED, CREATING DANGEROUS SITUATIONS FOR PEOPLE OR THINGS.

WHEN IN OPERATION, THE DEVICE MUST BE LOCATED AT A SAFE DISTANCE FROM ALL FLAMMABLE OBJECTS.

AT THE END OF USE IT IS NECESSARY TO TURN OFF THE DEVICE AND DISCONNECT IT FROM THE ELECTRICITY NETWORK. NEVER LEAVE THE DEVICE USELESSLY CONNECTED, OR TURNED ON WITHOUT SURVEILLANCE.

This product can be used by children over 8 years of age and people with reduced physical, mental or sensory abilities, or lack of experience and knowledge only under supervision or if they have received instructions regarding the safe use of the equipment and understand the risks associated with its use.

- It is NOT a toy, keep out of the reach of children and people with reduced psycho-physical abilities. DO NOT leave the product unattended in the presence of children or people with reduced mental and physical abilities. Do not allow children to play near the product.
- Always use on a firm, level, clear surface, making sure all parts are in contact with the floor. Keep away from streams, heights or railings. Do not use sharp objects near the product.
- Choose an adequate space where to place the device in order to guarantee free access and the best safety conditions. Leave a free space of at least 100 cm on both sides of the equipment and a minimum space of 100 cm in front.
- Before each exercise session, make sure there are no foreign objects under the vibration plate, and remove any objects promptly. Do not use the tool if there are objects under it.
- Make sure that the environment is adequately ventilated and with optimal oxygen exchange during use. Avoid drafts.
- Do not expose to bad weather, heat sources or open flames. Avoid humid environments or near water. This device must remain dry. Avoid prolonged exposure to direct sunlight. Do not make any structural changes to the product.
- Use the product one person at a time. Before each workout, always stretch. If you experience discomfort, chest pain, or any other symptoms, discontinue use and consult a physician. If you feel tired, discontinue use and rest.
- The use of this device is not for therapeutic purposes.
- A weighted use of this tool is recommended, in compliance with one's psychophysical conditions, away from meals and from the intake of medicines or substances that can in any way affect balance, vision, stability and in general compromise the general physical state. Do not use during pregnancy.
- Do not use the platform in case of previous, temporary, suspected or confirmed pathologies, problems related to balance, pregnancy, use of drugs, alcohol and medicines.
- Do not use in case of thrombosis, phlebitis or other pathologies affecting the venous, arterial and/or lymphatic system. Do not use in case of swelling or acute inflammation. Do not use in case of heart, cardiovascular or orthopedic problems.


- The platform is not designed for use by children. Children under the age of 12 must not use the tool. Young people aged 12 and over can use the tool under the constant supervision of an adult.
- Always wear appropriate clothing; do not wear clothing that could in any way hinder use. Always make sure that hands, hair and clothing are not trapped by the moving parts of the appliance. Do not use with bare feet, always wear only footwear suitable for the activity. Do not use the product if you are in an impaired state.
- Maximum user weight: 100kg. Don't overload. Do not use for purposes other than those intended. Avoid using in case of breakage or defects.
- Don't move when I load. Disconnect from the power supply before any movement. Do not leave at the disposal of children or pets.
- To use this tool, you must be in perfect psychophysical condition. A check-up with your doctor is recommended before using the device.
- The remote works with 2 AAA batteries (NOT included in the package). In the event of a battery fluid leak, do not allow the leaked fluid to contact your skin or clothes. If battery fluid gets on your skin or clothing, wash the skin with water immediately. **ATTENTION!** Danger of explosion if battery is replaced by an incorrect type. Only use and replace batteries of the same size and type. Improper use of the battery may cause leakage, overheating or explosion. This liquid is corrosive and can be toxic. It can cause skin and eye burns and is harmful to ingest. **DO NOT** heat, open, puncture, cut or throw the battery into a fire. Batteries are special waste; therefore, they cannot be mixed with common rubbish but must be disposed of separately, in the appropriate collection containers, according to current regulations. Before inserting the new battery, thoroughly clean the battery compartment with a soft, dry cloth, or follow the battery manufacturer's recommendations for cleaning. Do not immerse in water, protect from splashes, keep away from heat sources, direct sunlight or open flames. Do not reverse the polarity of the batteries. Remove the batteries if the device is not to be used for an extended period. Do not ingest, keep out of reach of children.

ASSEMBLY AND USE:

The device comes already assembled.

Positioning of the vibration plate:

1. Place the device on a solid, flat and clear surface, making sure that all parts are in contact with the floor. Keep away from streams, heights or railings. Do not use sharp objects near the product.
2. Choose an adequate space where to place the device in order to guarantee free access and the best safety conditions. Leave a free space of at least 2m from the ceiling, 1m on both sides of the apparatus and a minimum space of 1m in front. The free area around the device must include the space necessary for a possible emergency descent.
3. Before each exercise session, make sure there are no foreign objects under the vibration plate, and remove any objects promptly. Do not use the tool if there are objects under it

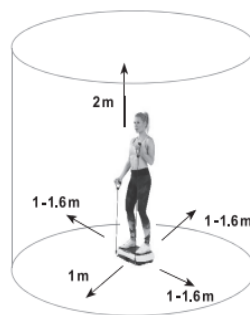Connection of the accessories:

If you also want to exercise the upper part of the body intensely (e.g. biceps, triceps, etc.) you can use the elastic bands supplied. For this purpose, it is sufficient to hook the carabiners present at the ends of the elastic bands to the specific hooks positioned at the base of the platform.


ATTENTION: if you do not want to use the elastic bands, it is recommended to keep them unhooked from their housings and to keep them in a safe place

INSTRUCTIONS FOR USE:

ATTENTION: When using for the first time, it is recommended to start training from the minimum vibration level. Improper use of the device can be dangerous to your health. Do not jump, kneel or stamp your feet on the device when it is switched on or off. Do not stand only on one side of the platform, but place your feet parallel, creating a stable and secure footing. The device can be operated using the buttons on the control panel or the remote control.

To turn on the device:

- Always plug the power cord into the tool before plugging it into the power outlet. After connecting the cable in this order, turn on the tool by pressing the main switch (4).
- Before disconnecting the cable from the socket, switch off the tool by pressing the main switch. Then disconnect the cable from the socket, and only then remove the power cable from the tool. It is recommended to turn off the main switch before removing the power cord.
- Once the device has been turned on, the displays (8) and (10) will show “HH” and “00” respectively. The platform will be in manual mode: the duration of the training and the level of vibration can be selected independently. To select the duration of the workout, press the TIME+ (7 or 19) and TIME- (17 or 27) buttons. The maximum training time is 15m. To activate the vibration, press the START button (16) on the platform or ► (18) on the remote control. After a 3 sec countdown, the treadmill will begin to vibrate. The vibration speed can be increased or decreased using the SPEED+ (11 or 28) and SPEED- (12 or 20) buttons.
- The platform provides three vibration zones depending on the training you wish to undertake: DELICATE (central zone), MEDIUM (middle zone), INTENSIVE (external zone). Place your feet on the platform according to the vibration you want to receive:

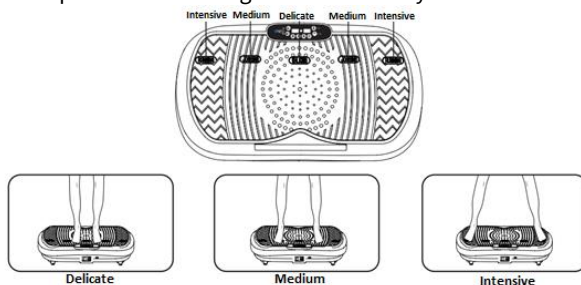- The device will automatically turn off when the workout time is over. The activity can also be interrupted at any time by pressing the STOP key (15) on the platform or (25) on the remote control.
- Manual mode can also be selected using the MANU key (14) on the control panel or M (21) on the remote control.
- To train with a predefined program (from P0 to P9), choose one of the ten programs by repeatedly pressing the AUTO button (13 or 26) until the desired program is selected. All programs last 15 minutes and the vibration level is automatically changed every minute by the device. When a predefined program is selected, the displays will show one the current vibration level, the other the program number and the minute countdown alternately. The 10 predefined programs are described in detail in the following table:

Program	MIN 15	MIN 14	MIN 13	MIN 12	MIN 11	MIN 10	MIN 9	MIN 8	MIN 7	MIN 6	MIN 5	MIN 4	MIN 3	MIN 2	MIN 1
P0	04	08	10	12	17	20	24	20	17	06	12	17	20	24	20
P1	02	04	06	08	10	13	15	12	09	08	08	10	13	15	12
P2	08	14	08	14	08	20	08	14	08	14	14	08	20	08	14
P3	12	23	29	38	49	56	68	76	48	29	38	49	56	68	76
P4	09	13	16	21	24	35	67	50	38	29	21	24	35	67	60
P5	23	55	23	55	23	55	23	55	23	55	55	23	35	23	55
P6	11	11	25	25	38	38	48	48	25	25	25	38	38	48	48
P7	30	30	50	50	70	70	40	40	60	60	50	70	70	40	40
P8	55	40	55	40	50	77	80	65	65	99	40	50	77	80	65
P9	40	85	75	65	25	96	78	60	65	80	65	25	96	78	60

Programs:

- Standby mode can be activated by pressing the relevant button (29) on the remote control. In standby mode, the displays will automatically turn off.

TRAINING TIPS:

Upright position: useful for loosening up muscles and stimulating circulation.

Push-ups: for an intensive training of the shoulder and arm muscles.

Seated stretch: Helps to relax the upper body and hips and work out the shoulder muscles and abs.

Leg relaxer: useful for loosening and relaxing the muscles of the legs and hips.

Training with elastic bands:

Lateral pull-ups: in an upright position, grasp the handles of the bands and raise both arms until they are parallel to the ground. Then lower your arms to the starting position

Crossed tractions: in an upright position, grasp the handles of the bands and lift one arm at a time until it is parallel to the ground, crossing the movement on the opposite side of the body and rotating the head in the direction of the movement.

Back Step/Bicep Curl: Standing with one foot on the platform, grasp the handles of the bands and bend the knee of the leg not resting on the platform. Keeping your arms close to

your torso, bend your elbows as you lift your hands towards your chin. Stand up by lowering your arms.

Push-up/front pull-up: From a standing position, grasp the handles of the bands and perform a push-up. At the same time raise your arms keeping them straight forward. Raise yourself to a standing position by lowering your arms.

Standing Biceps Curl: While standing, grasp the handles of the bands and bend both arms until your hands reach your shoulders. Then lower your arms to the starting position.

Alternatively, it is possible to perform the same exercises without the help of the training bands, to be replaced with dumbbells, bottles of water, or other weights (accessories not included in the package).

USE

Warnings

If any unexpected anomaly occurs during the use phase of the device, it is recommended to follow the procedure below:

- Turn it off.
- Disconnect it from the mains.
- Contact the manufacturer immediately.

In case of any anomaly of any electromagnetic nature and not foreseen after having powered the device for domestic use, it is recommended to follow the procedure below:

- Turn it off.
- Disconnect it from the mains.
- Check that the configuration has been carried out correctly.
- If there are no anomalies in the configuration, but the problem persists, contact the manufacturer immediately.

This device is not a toy, therefore it must never be made available to children.

It is forbidden to use the device during the hours of rest. It is not permitted to keep the device switched on continuously. It is recommended to take breaks to avoid overheating. Never cover the device, especially when it is turned on. If it stops working, malfunctions or makes unusual noises, turn it off using the switch and disconnect the plug from the electrical socket. Consult a qualified technician or the manufacturer. It is absolutely forbidden to force the operating mechanism of the product in any shape or form.

Periodic maintenance

Cleaning and maintenance operations must not be carried out by children unless they are older than 8 and are supervised. Keep the appliance and its cord out of the reach of children under 8 years of age.

Periodically it is necessary to clean the device from any accumulation of dust or other dirt that may have formed and/or deposited on its external surface. If necessary, use a compressor working with a commensurate pressure so as not to damage its components. Then, or in the absence of the compressor, use a brush with very soft bristles to avoid scratching or damaging the surfaces. It is also possible to use a damp cloth possibly accompanied by a multipurpose, non-abrasive product, free of alcohol or aggressive

solvents. Do not use abrasive sponges and chemical solvents. During cleaning, prevent water from touching internal electrical parts.

STORAGE

During periods of non-use, the device must be stored in a dry and clean place, possibly in the original packaging.

PUTTING OUT OF SERVICE

The device is manufactured and built according to criteria of robustness, durability and flexibility that allow it to be used for a relatively long period. This period may vary based on conditions and frequency of use.

This device is classified as W.E.E.E. (electrical and electronic equipment), once it reaches the end of its technical and operational life it must be taken out of service in conditions that it can no longer be used for the purposes for which it was designed and built at the time, making it possible to reuse it. of the raw materials that constitute it. Contact the relevant collection center for instructions on decommissioning the device.

The same deactivation procedures must be observed in all the following cases:

- Decommissioning of the device and storage in warehouse.
- Final dismantling and subsequent disposal.

THE MANUFACTURING COMPANY DOES NOT ASSUME ANY LIABILITY FOR DAMAGES TO PERSONS OR THINGS RESULTING FROM THE REUSE OF SINGLE PARTS OF THE DEVICE FOR FUNCTIONS OR ASSEMBLY CONFIGURATIONS DIFFERENT FROM THE ORIGINAL THOSE. THE MANUFACTURING COMPANY REJECTS ANY ACKNOWLEDGMENT, IMPLIED OR EXPRESS, OF FITNESS FOR A SPECIFIC PURPOSE OF PARTS OF THE DEVICE REUSED AFTER DEACTIVATION FOR ITS DISPOSAL

DISCARDING WORN HOUSEHOLD DEVICES

The symbol indicated in Legislative Decree no.49 of 14 March 2014 has been reported on the

product label, which transposes directive 2015/863 / EU, which indicates the need to dispose of the product and its components in separate waste, specifically in electrical and electronic waste. Contact the waste collection center in your area for technical support

THE MANUFACTURER IS IN NO WAY RESPONSIBLE FOR DAMAGE CAUSED BY THE DEVICE IF NOT USED IN THE FULL VERSION AND FOR THE USES AND METHOD OF USE SPECIFIED IN THIS MANUAL. THE MANUFACTURER IS IN NO WAY RESPONSIBLE FOR ANY DAMAGE TO PERSONS OR THINGS ARISING FROM THE RECOVERY OF USED PARTS OF THE DEVICE AFTER ITS DISPOSAL.

The paper and cardboard used in the appliance packaging are recyclable.

Under no circumstances should the packaging be disposed of in common rubbish or disposed of in the environment. Packaging must be disposed of in a differentiated way, separating them by type as per the regulations in force in your country, to facilitate the raw material recycling process.

Contact your municipality or an authorized waste collection center for relative and more detailed information.

For this appliance, the regulations relating to the recycling of consumer goods packaging by means of the National Packaging Consortium have been applied.

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The object of this declaration (same as point 3) is comply with harmonized EU regulation

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Description:

Vibrating platform with accessories and remote controller. With 10 programs and 99 levels of vibration.

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